



NITINOL STONECATCHER



NiTi 4xDØ - KØ - Form - x

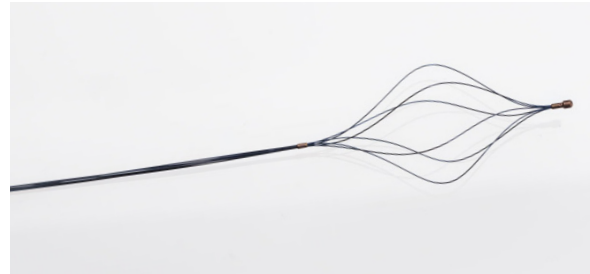

NiTi 4 (6) - DØ - KØ - Form

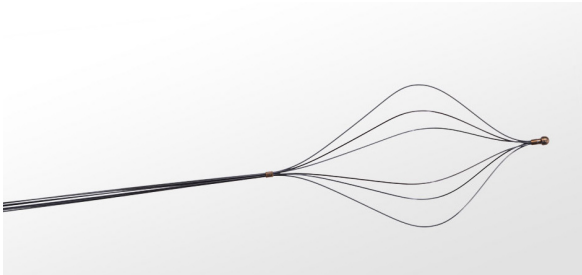

NiTi 4 (6) - DØ - KØ - Form

Explanations

DØ	diameter of wire
KØ	diameter of basket (15, 20, 25, 30, 35 mm)
Form	GF: twisted form DF: diamond form
x	Litho: lithotripsy capable Fü: guide wire capable (only KØ 20, 25, 30, 35)

Construction

Basket 4 or 6 Nitinol wires combined into a stainless steel tip.

Pull rope made of a stainless steel strand or continuous Nitinol wire (for lithotripsy capable basket) and stainless steel pull rod.

Customer requirements

In addition to our standards we also manufacture stonecatchers according to your specifications and wishes.

